

ARRAKIS Massage Chair User Manual



Contents

Contents	1
Safety Precautions	2
Usage Environment	2
Medical Warning	2
Preparation Before Use	3
1.How to install the cushion	3
2.How to adjust the position	3
3.How to move the chair	4
4.Massage chair placement	4
5.Power on	5
Main Structure and Function	6
Hand Controller Introduction	7
Program Operating Instruction	8-9
Daily Maintenance	10
Product Specification	10

Please read the user manual carefully before using the chair.

Daily Maintenance

- 1.Wipe off dust and dirty from the massage chair with a dry, soft, clean cloth after use.
- 2.Do not use water or other cleaning detergents to clean the massage chair.
- 3.Keep the massage chair clean and away from places of high humidity.
- 4.If the massage is not to be used for a long period of time, cover the massage chair with a dust cover and unplug its power plug from the electrical outlet.
- 5.The recommended duration of the operation is 30 minutes. It is not advisable for user to exceed the suggested time frame

Product Parameter

Model: **ARRAKIS Massage Chair**

Rated Voltage:110~250V,50/60Hz

Rated Power:90W

Rated Working Hours:15min

Test Description:Qualified

Certificate: CE

Remote Control Usage Guide



• Back Down

Press this button, the backrest will be moved down to the lowest position, release this button, it will stop.





• Massage Intensity

Under manual mode, Press “” button to increase the massage speed, Press “” button to decrease the massage speed



• Airbag Intensity

Open airbag function, Press “” button to increase the airbag intensity, Press “” button to decrease the airbag intensity



• Timing

Press this button to set the timing of 15-20-25-30 minute



• Zero Gravity

Press this button, the chair will go to zero-gravity position, Press it

Mode Switching

Auto Mode: In auto mode, there are 20 auto programs for choice.

Manual Mode: In manual mode, Choose any of manual function to have a specific massage.

Bluetooth Connection

Bluetooth function is turned on when the machine turned on

Turn on your mobile Bluetooth function search the machine Bluetooth and connection to playing the music.

Bluetooth Name: LEK-AMY Pin Code: 8888

Safety Precautions

Usage Environment

- To prevent the mildew, electric shock, fire or damaged, keep the chair away from wet place, Such as: Bathroom, Sauna, Hot spring area
- Its prohibited to use in an environment with room temperature above 40 degrees
- Do not use the chair in places where there is too much dust or lack of oxygen
- Do not put the chair in the hot stove or direct sunlight place, it may cause the accessories in the massage chair to overheat and cause fire
- Put in a clean place, nonsupport to use in dust, greasy dirt or oxygen lack area
- Put the chair in a flat place to use
- Make sure the socket is connected with the ground wire to prevent electric shock
- Do not use the chair when you are smoking


Medical Warnings

For following person please consult a doctor before using

- Malignant tumor patients
- People who have heart attacks
- Severe dermatosis patients
- Severe bone hyperplasia patients
- Pregnancy
- People who is in therapy under the instructions of a doctor
- The accident cause spinal abnormality or back bending people
- Joint deformation disease
- People was fitted with a pacemaker
- Visceral diseases people. such as: gastritis, esoenteritis, hepatitis. etc
- Severe blood circulation disorders
- High fever person
- Perception disorder patients
- Please take care of the children and do not let them play on the massage chair
- This product is not intended to be used by persons (including children) who are Limb insufficiency, sensory or mentally impaired, or lacking relevant experience and knowledge, unless someone responsible

Preparation Before Use

1.How to install the cushion

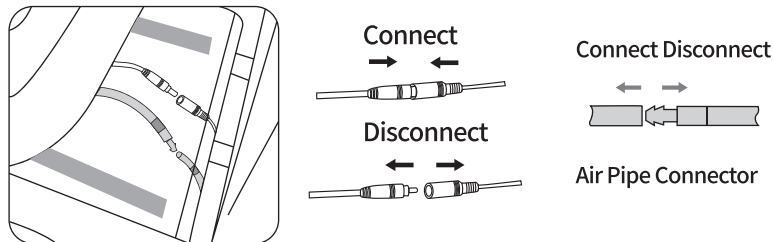

2.How to adjust the position














Press the Up & Down button to adjust the body to the most comfortable position

Please make sure there have enough space in front and behind the chair when you adjust the backrest and foot part



Program Operating Instructions

Remote Control Usage Guide

-  • On/Off
Press this button to start/stop the massage system.
-  • Auto Mode
Press this button to switch an auto program each time, there are 20 auto programs for choice
-  • Neck Massage
Press this button to start/stop the neck massage
-  • Back Massage
Press this button to start/stop the back massage
-  • Neck Direction
Press this button to choose the neck massage direction
-  • Back Direction
Press this button to choose the back massage direction
-  • Roller Speed
Press this button to switch foot roller speed each time ,3 levels of massage speed for choice
-  • Heating Function
Press this button to start/stop the heating in the back,seat and calf
-  • Foot Roller Massage
Press this button to start/stop the foot roller massage
-  • Head/Shoulder Airbag Massage
Press this button to start/stop head and shoulder airbag massage
-  • Arm/Seat Airbag Massage
Press this button to start/stop arm and seat airbag massage
-  • Calf Airbag Massage
Press this button to start/stop calf airbag massage
-  • Back Up
Press this button, the backrest will be moved up to the original position,release this button, it will stop.

Hand Controller Introduction



3.How to move the chair

1. Press down the back of massage chair to approximately 45° degree angle.

2. This will allow it to easily be moved to another location.

CAUTION:


* Do not move massage chair when occupied or in use.

* Moving the massage chair may cause damage to floors. Floor protection is recommended.

4.How to place the chair


Please leave enough space for in front and behind the chair

- The place where the massage chair is placed needs enough space for the reclining of the seat and the rise of the foot part
- Make sure there is at least 40cm behind the massage chair and at least 40cm in front of the chair
- Make sure there are no obstructions in the front and back of the massage chair
- Place a carpet under the massage chair to prevent the ground being damaged

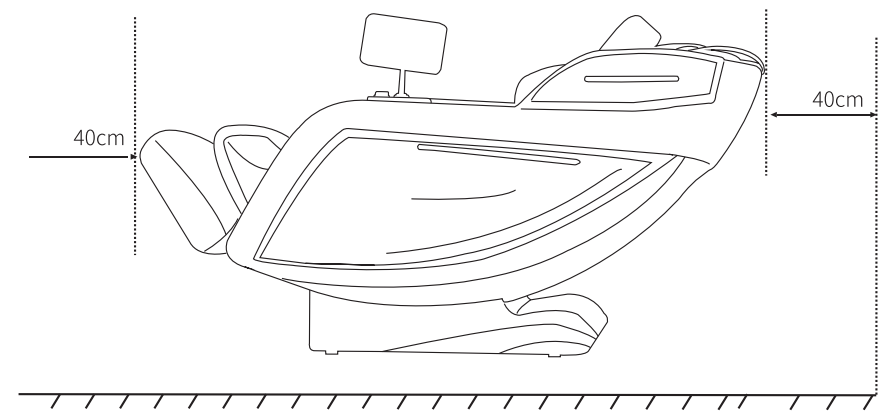

5.Power On

1.Make sure there are no obstacles around the massage chair

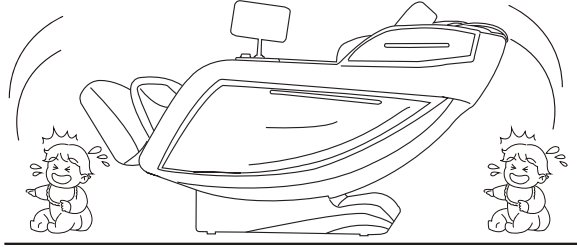


2.Insert the plug into the power outlet and turn on the power switch

*Check the power switch before plugging the massage chair into the outlet.

*The power switch should sit in the OFF position any time you unplug or plug in the chair.

*Follow all instructions for proper ground wire connectivity.

*Fully insert the plug into the outlet before turning the chair on.

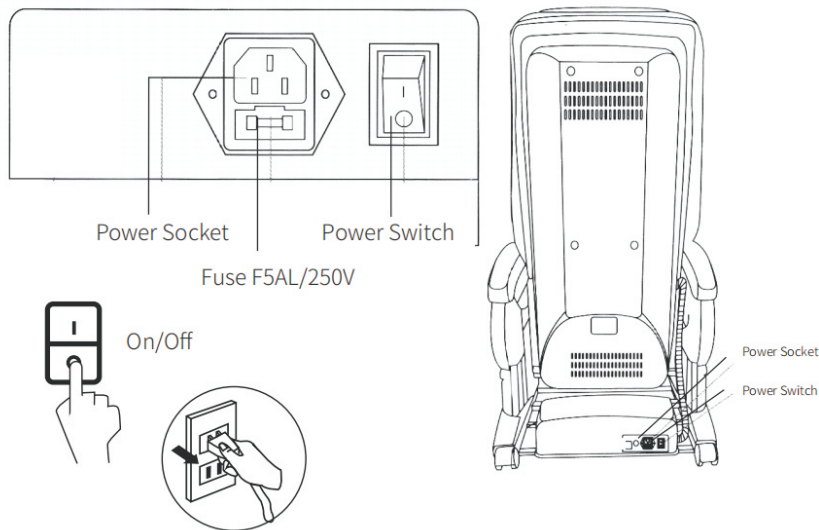

Main structure and function

